

Warning	1
Parts Introduction	2
Accessories	2
Basal Operation	3
40K Cavitation Explode Procedure	3
Caution for using 40K	4
Bipolar RF for facial	4
Tripolar RF for body	5
Hexpolar RF for body	5
Caution for using RF	5
Others	6
Forbidden group	6
Technical parameters	7
Accessories list	7

---

## Warning

Dear user:

Thanks for using our company's products, in order to fully make use of the products, we sincerely suggests follow:

- ◆◆ Please read the instruction in detail and keep it well for further reading and consulting.
- ◆◆ Please follow the instruction to fit and operate the machine correctly.
- ◆◆ Please do not remove and change any accessories of the machine.
- ◆◆ It forbids others to open and remove the mainframe except the authorized people in the company.
- ◆◆ It forbids others to change and maintain the machine except the professional people.
- ◆◆ Please contact us if there is any trouble during operation.

**This product does not belong to medical, please don't use for disease treatment!**

## Parts Introduction




- 1.LCD screen
- 2.IC card socket
- 3.General power switch
- 4.Power supply line jack
- 5.Hexpolar RF laser work head jack
- 6.40K cavitation explode work head jack
- 7.Tripolar RF work head jack
- 8.Bipolar RF work head jack
- 9.Service brake jack

## Accessories



Bipolar RF work head

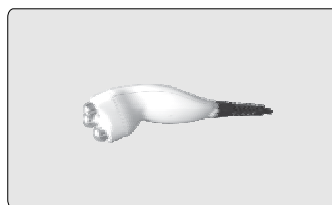

Tripolar RF work head



Hexapolar RF laser work head



40K cavitation explode work head



Service brake

## Basal Operation

1.Put the accessories well and connect them to the host.

- ①Insert Hexpolar RF laser work head to [5];
- ②Insert 40K cavitation explode work head to [6];
- ③Insert Tripolar RF work head to [7];
- ④Insert Bipolar RF work head to [8];
- ⑤Insert Service brake to [9].

2.Insert power supply line to [4], turn on the General power

switch [3], enter interface as (Fig 1) shows,input default Password "000000" to enter interface as (Fig 2).

3.Click "ENTER..." or click any place of the screen, enter Function interface as (Fig 3).



(Fig1)

## 40K Cavitation Explode Procedure

Apply the cavitation theory, uses the 40KHz ultrasound to make intense vibration, which causes a inner explosion in the fat tissues, thus get rid of lipocyte. It can effectively break down stubborn and deep fat, decrease lipocyte quantity and its size.

The diathermal of the low frequency ultrasound can promote lipocyte and tissues liquid flow, which help drainage of fatty acid, waste and toxin through lymph system.



1.Clean the parts and applied gel on fleshy parts

(e.g. abdomen, hip, thigh) .The amount depends on moving flexibility of the cavitation explode work head.

2.System default function is cavitation explode. Click the clock as (Fig 3) shows,<< time >> will be displayed and enter interface as (Fig 4), setting work time:

- (1)Click << and >> besides of the << time >> to set work time, press one time to change one minute;
- (2)Click 30:00 to change ten minutes per time;
- (3)Regulat work time range is 1:00~60:00 minutes, the system default time is 30 minutes.
- (4)Each part is ten minutes of every treatment.

3.Click << and >> besides of the << intensity >> to adjust output intensity: 5W(weak)~25W(strength).

4.Click  , the machine is in standby state. Beautician put the service brake under your feet, step on to start working, move away for pause. Click  , the machine stop working.

5.Beautician hold the cavitation explode work head to move slowly on skin, in circle or beeline repeatedly, and the Other hand push the fat towards the work head.



(Fig2)



(Fig3)

6. Do not use on back, avoid using on bones.
7. Do not use on uterus part when women do abdomen care. Abdomen treatment during menstruate should be avoided.
8. The machine will pause automatically if the setting time is over. Please click any place of the screen According to the hint.
9. Clean the leftovers with hot towel, use warm towel to clean the cavitation explode work head.

## Caution for using 40K

1. Applied enough gel on skin need treated.
2. Do not stay unmoved on any part of your body, do not stay on bones.
3. The total treatment time for daily should be within 30 minutes.
4. There will be ultrasonic only when you step on the service brake while using 40K cavitation explode work head.
5. Do not allowed to be placed under the state of no-load to avoid the damage caused by high heat.
6. 40K work head can be used on fresh part of the body, Do not use for facial. If the skin is sensitive or becomes unhealthy, the strength of the power should be adjusted down.

## Bipolar RF for facial

1. Click " Bipolar RF " as (Fig 3) shows, enter interface as (Fig 4).
  2. Set work time, 20 minutes for facial and 15 minutes for eyes treatment.
  3. Click << and >> besides of the << intensity >> to adjust Output intensity: 1(weak)~5(strength).
  4. Click (⏏), the machine is in standby state. Beautician put the service brake under your feet, step on to start working, move away for pause. Press (⏏) to stop working.
  5. Apply gel on face around skin, move the electrode probe in circle, from inside to outside, operating on both side, right and left, one side is over then to start the other side, and must reduce the intensity according to the less fat on forehead and eyes.
- NOTE :Do not let gel into eyes .
6. The machine will pause automatically if the setting time is over.
  7. Clean the leftovers with hot towel, use warm towel to Clean the work head.



(Fig4)

## Tripolar RF for body

1. Click “ Tripolar RF ” as (Fig 3) shows, enter interface as (Fig 5).
2. Set work time, 20 minutes for each part.
3. Click << and >> besides of the << intensity >> to adjust Output intensity: 1(weak)~10(strength).
4. Click ⏸, the machine is in standby state. Beautician put the service brake under your feet, step on to start working, move away for pause. Press ⏵ to stop working.
5. Applied gel on skin needed nursing, If there is too much fat, you can increase the output intensity. Tripolar RF work head suit for neck, arm, crus.
6. The machine will pause automatically if the setting time is over.
7. Clean the leftovers with hot towel, use warm towel to Clean the work head.



(Fig5)

## Hexapolar RF for body

1. Click “ Hexapolar RF ” as (Fig 3) shows, enter interface as (Fig 6).
  2. Set work time, 20 minutes for each part.
  3. Click << and >> besides of the << intensity >> to adjust Output intensity: 1(weak)~10(strength).
  4. Click ⏸, the machine is in standby state. Beautician put the service brake under your feet, step on to start working, move away for pause. Press ⏵ to stop working.
  5. The system default laser is start state, click 🔴, the laser function is closed.
- NOTE: Do not let the laser light straightly shine upon The eyes in case of eyes hurt.**
6. Applied gel on skin needed nursing, the electrode probe move slowly on skin, If there is too much fat, you can increase the output intensity.
  7. The machine will pause automatically if the setting time is over.
  8. Clean the leftovers with hot towel, use warm towel to Clean the work head.



(Fig6)

## Caution for using RF

- 1.Beautician should wear insulated gloves during operation.
- 2.Do not wear any aglet, both for the beautician and the users.
- 3.The electrode probe must be fully in contact with skin.
- 4.Do not make the RF electrode stay on any part without moving. Keep it moving on skin slowly.
- 5.Skin around eyes and forehead is very thin, the RF output intensity should be relatively weakened.
- 6.Do not contact skin with unipolar when using of bipolar, tripolar, hexpolar RF work head otherwise it will Produce pulse discharge, gel paste etc. Products like gel maybe evaporated into a mist.
- 7.Do not let the laser light straightly shine upon the eyes in case of eyes hurt.

## Others

Click “set” as the(Fig 2) shows. Enter interface as (Fig 7):

- ①Click“ basic set ”, enter interface as (Fig 8).
  - 1---Select language
  - 2---Regulate date and time
  - 3---Regulate brightness of the screen
  - 4---Click”calibration” enter interface as(Fig 9),for Calibrating the screen.
  - 5---Keys recovery back


**Notes “advanced setting” is for manufacturer only.**

## Forbidden group


- 1.People with heart disease or high blood pressure, or configured cardiac pacemaker.
- 2.Patients with acute inflammation, asthma, deep vein thrombosis, thyroncus, cancer.
- 3.People with hemorrhagic disease, trauma or who is bleeding.
- 4.Pregnant women
- 5.Medical Plastic parts, or parts with artificial in filling inside.
- 6.People with metal inside the body(include metallic tooth)
- 7.Patient with skin inflammation or with edema.
- 8.People with an abnormal immune system.
- 9.Numb or insensitive to heat.



(Fig7)

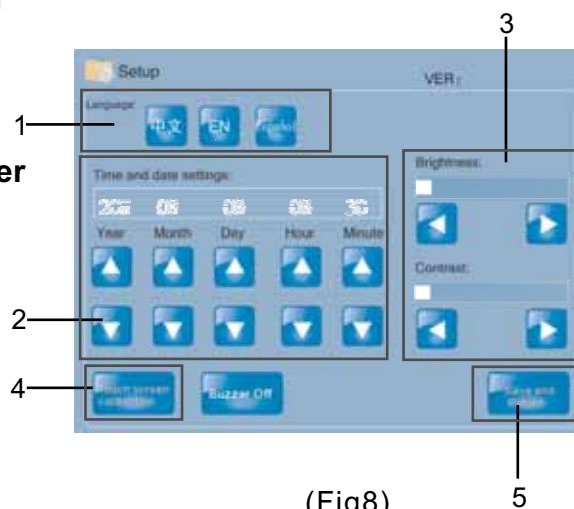

(Fig8)

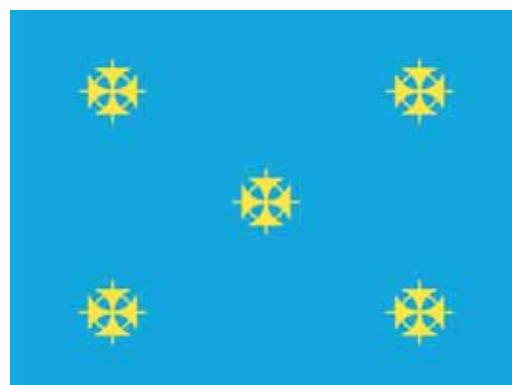

(Fig9)

---

## Technical parameters

Input: AC110V/60Hz ☐ AC220V/50Hz☐

Power:  $\leq 75\text{VA}$

Output RF: 5MHz

Output ultrasonic: 40KHz

## Accessories list

1.Host	1set
2.Power supply line	1pc
3.User manual	1copy
4.Bipolar RF work head	1pc
5.Tripolar RF work head	1pc
6.Hexpolar RF laser work head	1pc
7.40K cavitation explode work head	1pc
8.Service brake	1pc
9.Accessory frame	1pcs
10.Hexagonal screws	2pcs
11.Hexagon wrench	1pc